

UK Office

Everest Biotech Ltd Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries: info@everestbiotech.com Sales: sales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

US Office

Everest Biotech c/o Abcore 405 Maple Street, Suite A106 Ramona, CA 92065 USA

Inquiries: info@everestbiotech.com Sales: usasales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: 888-320-4628 (toll-free) Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB09003 - Goat Anti-AUTS2 Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: AUTS2, autism susceptibility candidate 2, KIAA0442, MGC13140, autism-related protein 1 Official Symbol: AUTS2 Accession Number(s): NP_056385.1; NP_001120703.1; Human GeneID(s): 26053 Non-Human GeneID(s): 319974 (mouse) Important Comments: This antibody is expected to recognize NP_056385.1 isoforms 1 and NP_001120703.1 isoform 2.

Immunogen

Peptide with sequence C-ERPPSHTLKDIEAR, from the C Terminus of the protein sequence according to NP_056385.1; NP_001120703.1;.

Please note the <u>peptide</u> is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested


Peptide ELISA: antibody detection limit dilution 1:32000.

Immunofluorescence: Strong expression of the protein seen in the nuclei and cytoplasm of A431 and U2OS cells. Recommended concentration: 10µg/ml.


Flow Cytometry: Flow cytometric analysis of A431 cells. Recommended concentration: 10ug/ml.

Species Reactivity


Tested: Human Expected from sequence similarity: Human, Mouse, Dog

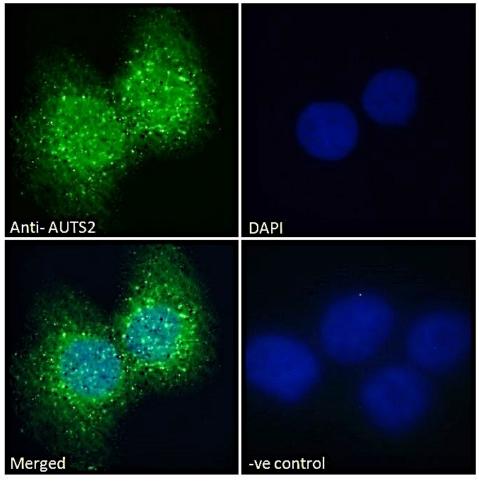

EB09003 Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

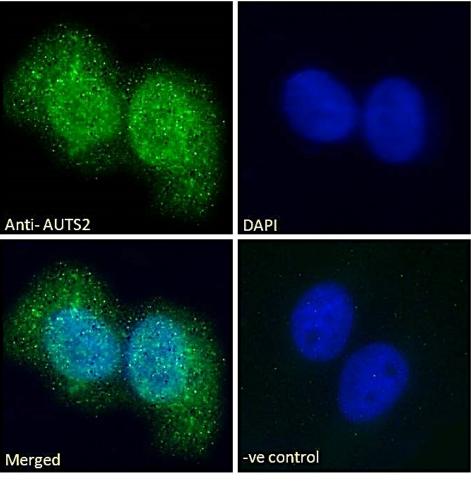

EB09003 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

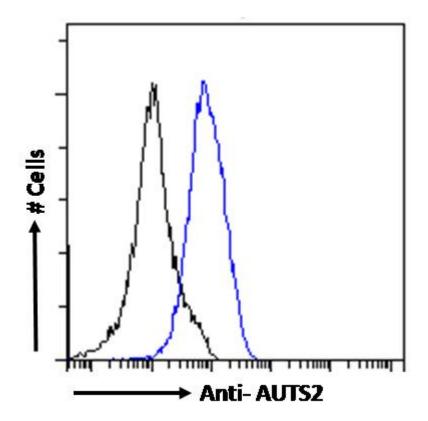

EB09003 Flow cytometric analysis of paraformaldehyde fixed A431 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control: Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.