

### **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre

77 Heyford Park
Upper Heyford
Oxfordshire
OX25 5HD

UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

### **US Office**

### **Everest Biotech c/o Abcore**

405 Maple Street, Suite A106

Ramona, CA 92065 USA

Inquiries:

info@everestbiotech.com

Sales:

 $\underline{usasales@everest biotech.com}$ 

Tech support:

support@everestbiotech.com

Tel: 888-320-4628 (toll-free)

Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB08191 - Goat Anti-ERCC1 Antibody

Size: 100µg specific antibody in 200µl



# **Target Protein**

**Principal Names:** ERCC1, excision repair cross-complementing rodent repair deficiency, complementation group 1 (includes overlapping antisense sequence), COFS4, UV20,

excision repair cross-complementing 1, excision repair protein

Official Symbol: ERCC1

Accession Number(s): NP\_973730.1; NP\_001974.1; NP\_001159521.1

Human GeneID(s): 2067

### **Immunogen**

Peptide with sequence C-QVDVKDPQQALKE, from the internal region of the protein sequence according to NP\_973730.1; NP\_001974.1; NP\_001159521.1.

Please note the peptide is available for sale.

### **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

# **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:4000.


Immunofluorescence: Strong expression of the protein seen in the nuclei of A431 and

U2OS cells. Recommended concentration:  $10\mu g/ml$ .


# **Species Reactivity**

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

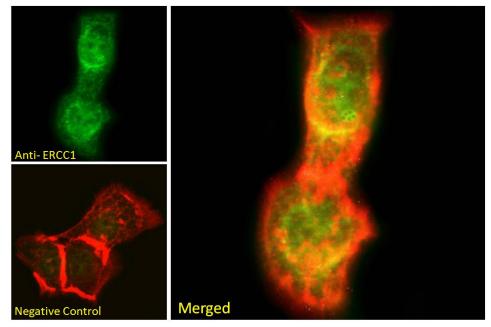

EB08191 Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. Actin filaments were stained with phalloidin (red). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

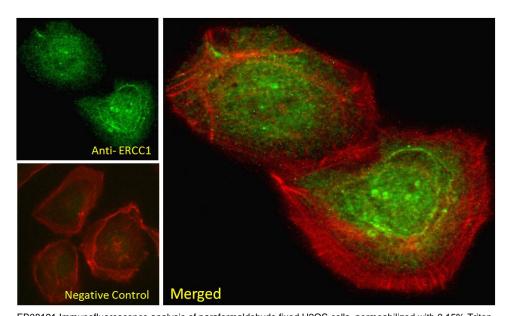

EB08191 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. Actin filaments were stained with phalloidin (red). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).