

### **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD

UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

#### **US Office**

#### **Everest Biotech c/o Abcore**

405 Maple Street, Suite A106 Ramona,

CA 92065 USA

Inquiries:

info@everestbiotech.com

Sales:

usasales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: 888-320-4628 (toll-free)

Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB06060 - Goat Anti-FBXO44 Antibody

Size: 100µg specific antibody in 200µl



### **Target Protein**

Principal Names: FBG3, FBX30, Fbx44, Fbxo6a, MGC14140, F-box gene 3, F-box

protein FBX30, F-box protein 44

Official Symbol: FBXO44

Accession Number(s): NP\_904319.1; NP\_001317284.1

Human GeneID(s): 93611

# **Immunogen**

Peptide with sequence C-AALTPPEPPSAEP, from the C Terminus of the protein sequence according to NP\_904319.1; NP\_001317284.1.

Please note the peptide is available for sale.

# **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum

Aliquot and store at -20°C. Minimize freezing and thawing.

# **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:8000.

**Western blot:** Preliminary testing showed a band at approx 26kDa in Human Kidney lysate after 2µg/ml antibody staining (calculated MW of 25.7kDa according to NP\_904319.1). Primary incubation 1 hour at room temperature.


**Immunofluorescence:** Strong expression of the protein seen in the nuclei of A431 and U2OS cells. Recommended concentration: 10µg/ml.

**Flow Cytometry:** Flow cytometric analysis of A431 cells. Recommended concentration: 10ug/ml.

## **Species Reactivity**


Tested: Human

Expected from sequence similarity: Human

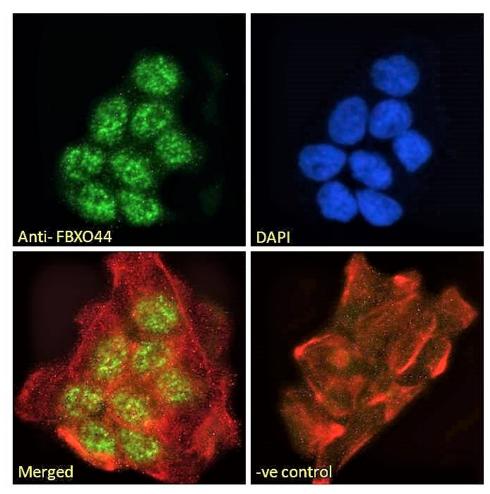

EB06060 Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. Actin filaments were stained with phalloidin (red) and the nuclear stain is DAPI (blue). Negative control:


Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

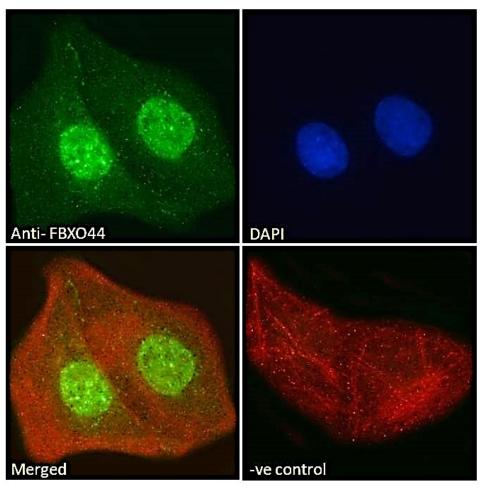

EB06060 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton.

Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. Actin filaments were stained with phalloidin (red) and the nuclear stain is DAPI (blue). Negative control:

Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

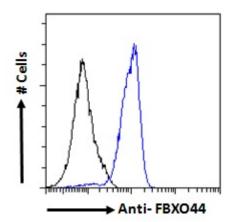

EB06060 Flow cytometric analysis of paraformaldehyde fixed A431 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.