

#### **UK Office**

### **Everest Biotech Ltd**

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD

UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB09143 - Goat Anti-SLC12A1 / NKCC Antibody

Size: 100µg specific antibody in 200µl



## **Target Protein**

**Principal Names:** SLC12A1, solute carrier family 12 (sodium/potassium/chloride transporters), member 1, BSC1, MGC48843, NKCC2, NKCC2A variant A, Na-K-2Cl

cotransporter, sodium potassium chloride cotransporter 2

Official Symbol: SLC12A1

Accession Number(s): NP\_000329.2

Human GeneID(s): 6557

## Immunogen

Peptide with sequence C-PHYEETSFGDEAQKR, from the internal region of the protein sequence according to NP\_000329.2.

Please note the peptide is available for sale.

## **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

## **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:32000.

**Western blot:** This anitbody has been successfully used in WB on Mouse: Alshahrani S et al. (2015) PMID: 26400961.

IHC: In paraffin embedded Human Kidney shows strong membrane staining in some CT cells. Recommended concentration, 3-6µg/ml.

**IF/ICC:** This antibody has been successfully used in IF on Mouse: Alshahrani S et al. (2015) PMID: 26400961. This antibody has been successfully used in ICC on Pig: Andreas Brännström et al. (2022) PMID: 36471055.

### **Species Reactivity**

Tested: Human, Mouse, Pig

Expected from sequence similarity: Human

## **Specific References**

### This antibody has been successfully used in ICC on Pig:

Andreas Brännström, Michael Hultström, Jenny Gustavsson, Zabih Aurfan & Mattias Günther

Intermittent thoracic resuscitative endovascular balloon occlusion of the aorta improves renal function compared to 60 min continuous application after porcine class III hemorrhage.

Eur J Trauma Emerg Surg (2022). https://doi.org/10.1007/s00068-022-02189-2 PMID: 36471055

#### This antibody has been successfully used in WB and IF in Mouse:


Alshahrani S, Almutairi MM, Kursan S, Dias-Junior E, Almiahuob MM, Aguilar-Bryan L, Di Fulvio M.

Increased Slc12a1 expression in  $\beta$ -cells

and improved glucose disposal in Slc12a2 heterozygous mice.

J Endocrinol. 2015 Dec;227(3):153-65.

PMID: 26400961

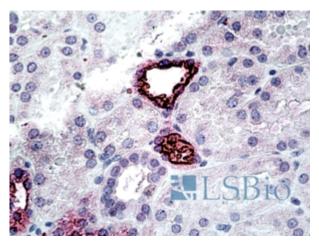

EB09143 (3.8ug/ml) staining of paraffin embedded Human Kidney. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.