

#### **UK Office**

**Everest Biotech Ltd** 

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD

UK

**Enquiries:** 

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB07921 - Goat Anti-TREM2 Antibody

Size: 100µg specific antibody in 200µl



**Target Protein** 

**Principal Names:** TREM2, triggering receptor expressed on myeloid cells 2, TREM-2, Trem2a, Trem2b, Trem2c, OTTHUMP00000039771, triggering receptor expressed on

myeloid cells 2a

Official Symbol: TREM2

Accession Number(s): NP\_061838.1

Human GenelD(s): 54209

#### **Immunogen**

Peptide with sequence C-HGQKPGTHPPSELD, from the C Terminus of the protein sequence according to NP\_061838.1.

Please note the <u>peptide</u> is available for sale.

## **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

### **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:128000.

IHC: In paraffin embedded Human Lung shows staining of alveolar macrophages.

Recommended concentration: 2-4µg/ml.

Immunofluorescence: Strong expression of the protein seen in U937 and HepG2 cells.

Recommended concentration: 10µg/ml.


**Flow Cytometry:** Flow cytometric analysis of U937 cells. Recommended concentration:

10ug/ml.


### **Species Reactivity**

Tested: Human


Expected from sequence similarity: Human

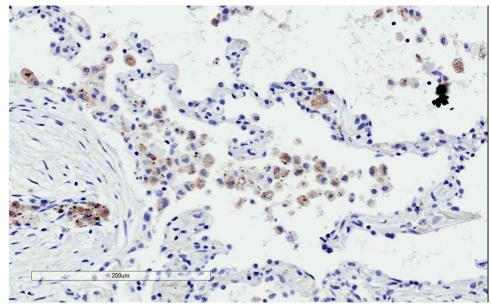


EB07921 (2µg/ml) staining of paraffin embedded Human Lung. Microwaved antigen retrieval with citrate buffer pH 6, HRP-staining.

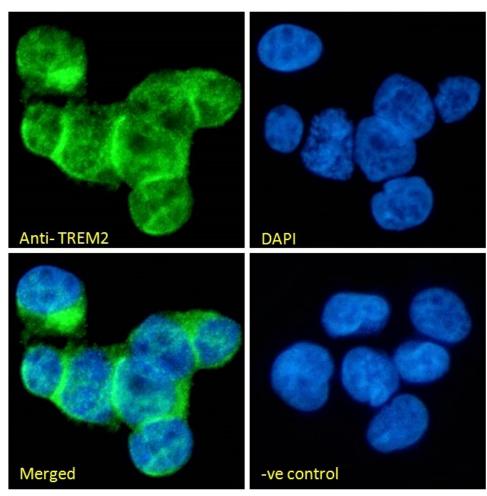

EB07921 Immunofluorescence analysis of paraformaldehyde fixed U937 cells immobilized on ShifixTM coverslip, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic and plasma membrane staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

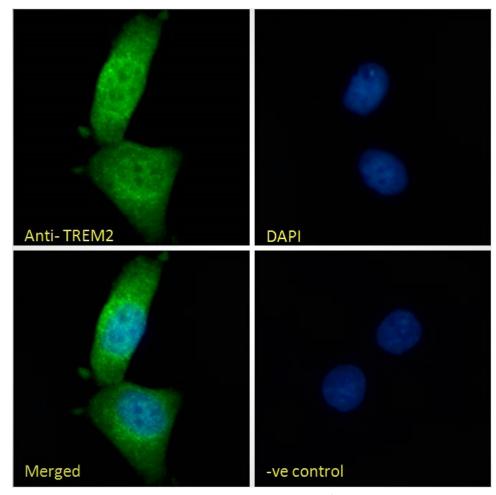

EB07921 Immunofluorescence analysis of paraformaldehyde fixed HepG2 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic and some nuclear staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

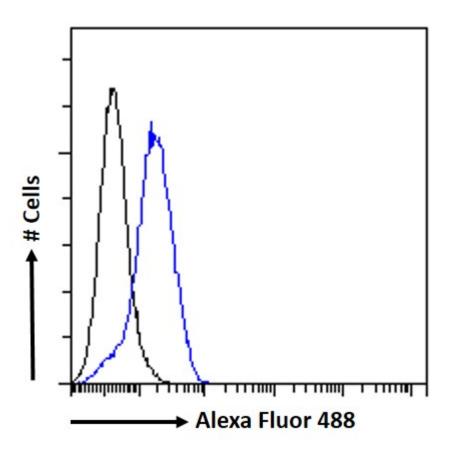

EB07921 Flow cytometric analysis of paraformaldehyde fixed U937 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.