

UK Office

Everest Biotech Ltd

Cherwell Innovation Centre
77 Heyford Park
Upper Heyford
Oxfordshire
OX25 5HD
UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

**Research Use Only. Not for
diagnostic or therapeutic use.**

EB07133 - Goat Anti-ADRB1 Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: ADRB1, adrenergic, beta-1-, receptor, ADRB1R, B1AR, BETA1AR, RHR, beta-1-adrenergic receptor

Official Symbol: ADRB1

Accession Number(s): NP_000675.1

Human GeneID(s): [153](#)

Immunogen

Peptide with sequence ESDEARRCYNDPK, from the internal region of the protein sequence according to NP_000675.1.

Please note the [peptide](#) is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:4000.

Western blot: This antibody is not recommended for use in Western blotting.

IHC: Paraffin embedded Human Heart. Recommended concentration: 4µg/ml.

Species Reactivity

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat, Dog

Specific References

This antibody (previous batch) has been successfully used in the following paper:

Ting-Yu Lin, Quynh N Mai, Hao Zhang, Emily Wilson, Huan-Chieh Chien, Sook Wah Yee, Kathleen M Giacomini, Jeffrey E Olgin, Roshanak Irannejad
Cardiac contraction and relaxation are regulated by distinct subcellular cAMP pools.
Nat Chem Biol. 2023 Jul 20.
PMID: 37474759

This antibody (previous batch) has been successfully used in Assay

Standardization of Rat Cardiomyocytes:

Gerd Wallukat, Harald Pruss, Johannes Muller, Ingolf Schimke
Functional autoantibodies in patients with different forms of dementia
PLoS One. 2018 Mar 14;13(3):e0192778.
PMID: 29538413

This antibody (previous batch) has been successfully used in the following paper:

Wenzel K, Schulze-Rothe S, Müller J, Wallukat G, Haberland A.
Difference between beta1-adrenoceptor autoantibodies of human and animal origin-Limitations detecting beta1-adrenoceptor autoantibodies using peptide based ELISA technology.
PLoS One. 2018 Feb 9;13(2):e0192615
PMID: 29425252

This antibody (previous batch) has been successfully used in the following paper:

Wenzel K, Schulze-Rothe S, Haberland A, Müller J, Wallukat G, Davideit H
Performance and in-house validation of a bioassay for the determination of beta1-auto
antibodies found in patients with cardiomyopathy.
Heliyon. 2017 Jul 31;3(7): e00362.
PMID: 28795160

This antibody (previous batch) has been successfully used in the following paper:

Annekathrin Haberland, Martin Holtzhauer, Alice Schlichtiger, Sabine Bartel, Ingolf
Schimke, Johannes Müller, Michael Dandel, Peter B Lippa, Gerd Wallukat
Aptamer BC 007 - A broad spectrum neutralizer of pathogenic autoantibodies against
G-protein-coupled receptors
Eur J Pharmacol. 2016 Oct 15;789:37-45.
PMID: 27375076

**This antibody (previous batch) has been successfully used for In Vitro Testing of
Apheresis Column Functionality in Rat:**


Wallukat G, Haberland A, Berg S, Schulz A, Freyse EJ, Dahmen C, Kage A, Dandel M,
Vetter R, Salzsieder E, Kreutz R, Schimke I.
The first aptamer-apheresis column specifically for clearing blood of β 1-receptor
autoantibodies.
Circ J. 2012;76(10):2449-55.
PMID: 22850243

**This antibody (previous batch) has been successfully used in IP on DNA aptamers
in vitro:**


Haberland A, Wallukat G, Dahmen C, Kage A, Schimke I.
Aptamer neutralization of beta1-adrenoceptor autoantibodies isolated from patients with
cardiomyopathies.
Circ Res. 2011 Oct 14;109(9):986-92.
PMID: 21868696

**This antibody (previous batch) has been successfully used in Assay
Standardization (beating frequency) of Rat Cardiomyocytes:**

Wallukat G, Muñoz Saravia SG, Haberland A, Bartel S, Araujo R, Valda G, Duchon D,
Diaz Ramirez I, Borges AC, Schimke I.
Distinct patterns of autoantibodies against G-protein-coupled receptors in Chagas'
cardiomyopathy and megacolon. Their potential impact for early risk assessment in
asymptomatic Chagas' patients.
J Am Coll Cardiol. 2010 Feb 2;55(5):463-8.
PMID: 20117461

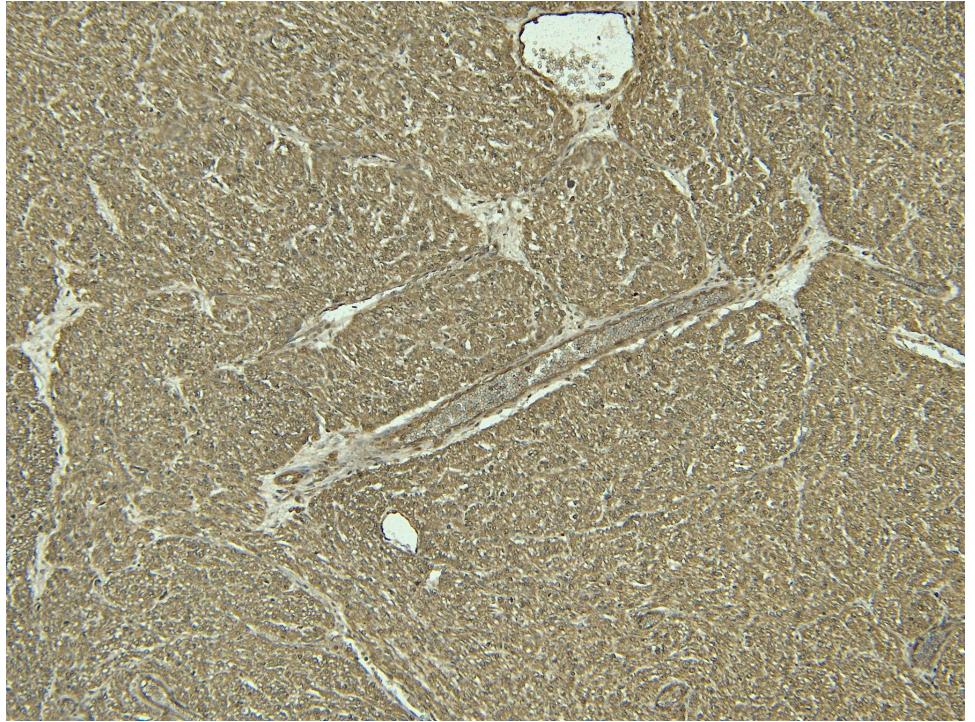

EB07133 (4 μ g/ml) staining of paraffin embedded Human Heart. Heat induced antigen retrieval with citrate buffer pH 6, HRP-staining.

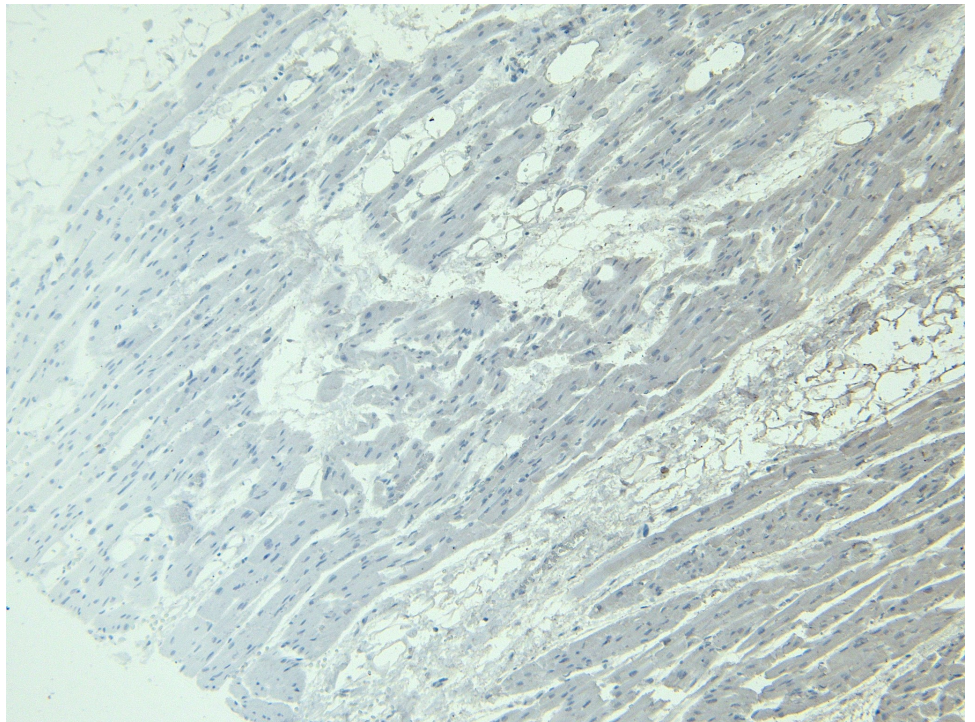

EB07133 Negative Control showing staining of paraffin embedded Human Heart, with no primary antibody.