

UK Office

Everest Biotech Ltd Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries: info@everestbiotech.com Sales: sales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

US Office

Everest Biotech c/o Abcore

405 Maple Street, Suite A106 Ramona, CA 92065 USA

Inquiries: info@everestbiotech.com Sales: usasales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: 888-320-4628 (toll-free) Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB06122 - Goat Anti-CENTB1 / ACAP1 Antibody

Size: 100µg specific antibody in 200µl

C

Target Protein

Principal Names: ACAP1, CENTB1, centaurin, beta 1, KIAA0050, KIAA0050 gene product, Arf GAP with coiled coil, ANK repeat and PH domains 1, ArfGAP with coiled-coil, ankyrin repeat and PH domains 1, centaurin, beta 1 Official Symbol: ACAP1 Accession Number(s): NP_055531.1 Human GeneID(s): <u>9744</u> Non-Human GeneID(s): 216859 (mouse)

Immunogen

Peptide with sequence C-EKLSRRSHDLHTL, from the C Terminus of the protein sequence according to NP_055531.1.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide. Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin. Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:32000. **Western blot:** Approx 80kDa band observed in Jurkat lysates (calculated MW of 81.5kDa according to NP_055531.1). Recommended for use at 0.3-1µg/ml.

Additional validation: This antibody has been successfully used in the following paper: Sikorski et al. (2018) PMID: 30377371.


Species Reactivity

Tested: Human Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

Specific Reference

This antibody has been successfully used in the following paper:

Krzysztof Sikorski, Adi Mehta, Marit Inngjerdingen, Flourina Thakor, Simon Kling, Tomas Kalina, Tuula A. Nyman, Maria Ekman Stensland, Wei Zhou, Gustavo A. De Souza, Lars Holden, Jan Stuchly, Markus Templin and Fridtjof Lund-Johansen A high-throughput pipeline for validation of antibodies Nat Methods. 2018 Nov;15(11):909-912 PMID: 30377371

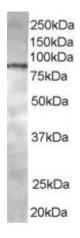

EB056122 staining (0.5µg/ml) of Jurkat lysate (RIPA buffer, 35µg total protein per lane). Primary incubated for 1 hour. Detected by western blot using chemiluminescence.