

International Office

Everest Biotech Ltd Vector Laboratories, Inc.

6737 Mowry Ave Newark, CA 94560 United States

Customer Service: <u>customerservice@vectorlabs.com</u> Technical Service: <u>technical@vectorlabs.com</u>

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB08605 - Goat Anti-CRIPTO Antibody

Size: 100µg specific antibody in 200µl

C

Target Protein

Principal Names: TDGF1, teratocarcinoma-derived growth factor 1, CR, CRGF, CRIPTO, Cripto-1 **Official Symbol:** TDGF1

Accession Number(s): NP_003203.1; NP_001167607.1 Human GenelD(s): <u>6997</u> Non-Human GenelD(s): 21667 (mouse), 680246 (rat)

Immunogen

Peptide with sequence YGRNCEHDVRKEN, from the internal region of the protein sequence according to NP_003203.1; NP_001167607.1.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested


Peptide ELISA: antibody detection limit dilution 1:4000.

Immunofluorescence: Strong expression of the protein seen in HepG2 and U2OS cells. Recommended concentration: 10µg/ml.


Flow Cytometry: Flow cytometric analysis of HepG2 cells. Recommended concentration: 10ug/ml.

Species Reactivity


Tested: Human Expected from sequence similarity: Human

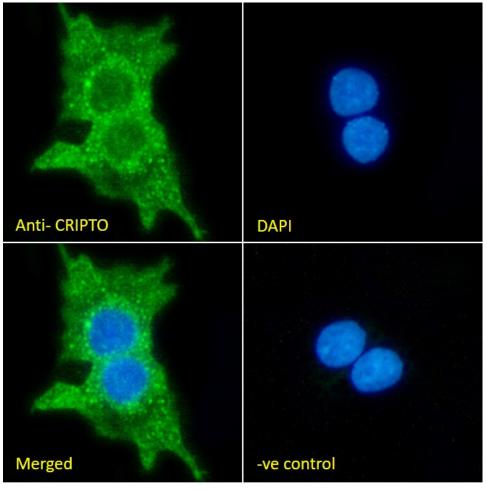

EB08605 Immunofluorescence analysis of paraformaldehyde fixed HepG2 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing vesicle and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

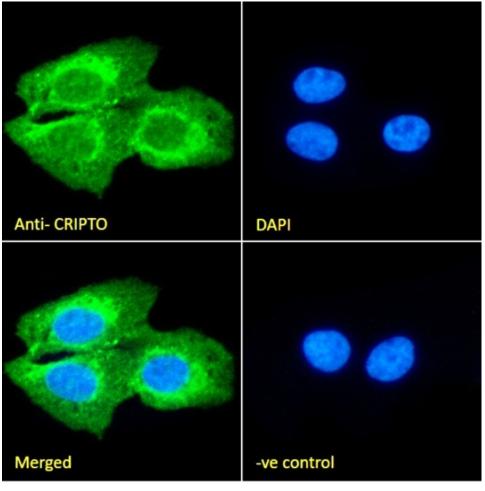

EB08605 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

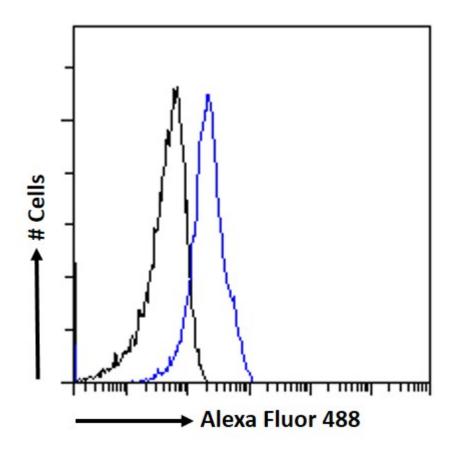

EB08605 Flow cytometric analysis of paraformaldehyde fixed HepG2 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control: Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.