

UK Office

Everest Biotech Ltd Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries: info@everestbiotech.com Sales: sales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

US Office

Everest Biotech c/o Abcore

405 Maple Street, Suite A106 Ramona, CA 92065 USA

Inquiries: info@everestbiotech.com Sales: usasales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: 888-320-4628 (toll-free) Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB09267 - Goat Anti-DAP12 Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: TYROBP, TYRO protein tyrosine kinase binding protein, DAP12, KARAP, PLOSL, DNAX-activation protein 12, killer activating receptor associated protein Official Symbol: TYROBP Accession Number(s): NP_003323.1; NP_937758.1

Human GenelD(s): 7305 Important Comments: This antibody is expected to recognize both reported isoforms

Immunogen

Peptide with sequence C-KQRITETESPYQE, from the internal region of the protein sequence according to NP_003323.1; NP_937758.1.

Please note the peptide is available for sale.

Purification and Storage

(NP_003323.1; NP_937758.1)


Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide. Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin. Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested


Peptide ELISA: antibody detection limit dilution 1:32000.IHC: Paraffin embedded Human Uterus and Spleen. Recommended concentration: 3-6µg/ml.

Species Reactivity

Tested: Human Expected from sequence similarity: Human

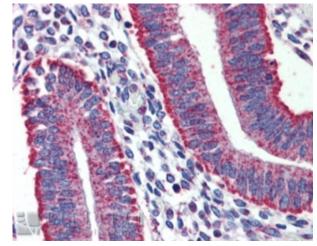

EB09267 (3µg/ml) staining of paraffin embedded Human Uterus. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.

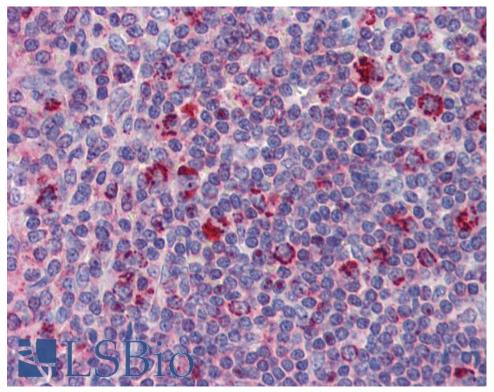

EB09267 (3µg/ml) staining of paraffin embedded Human Spleen. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.