

International Office

Everest Biotech Ltd

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

Customer Service:

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB05758 - Goat Anti-E2-EPF / UBE2S Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: E2-EPF, ubiquitin carrier protein, ubiquitin-protein ligase, ubiquitin-conjugating enzyme E2-24 kD, UBE2S, ubiquitin-conjugating enzyme E2S,

E2EPF, EPF5, ubiquitin carrier protein S, ubiquitin-protein ligase S

Official Symbol: UBE2S

Accession Number(s): NP_055316.2

Human GeneID(s): 27338

Non-Human GenelD(s): 77891 (mouse)

Immunogen

Peptide with sequence NSNVENLPPHIIR-C, from the N Terminus of the protein sequence according to NP_055316.2.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:32000.


Immunofluorescence: Strong expression of the protein seen in the cytoplasm and nuclei of A431 and HeLa cells. Recommended concentration: 10µg/ml.

Flow Cytometry: Flow cytometric analysis of HeLa cells. Recommended concentration: 10ug/ml.


Species Reactivity

Tested: Human


Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

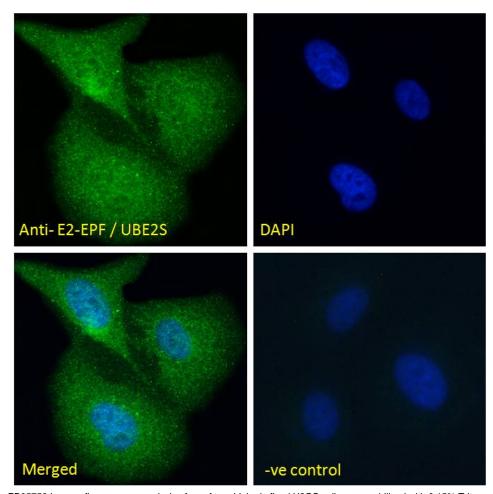

EB05758 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic and nuclear staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

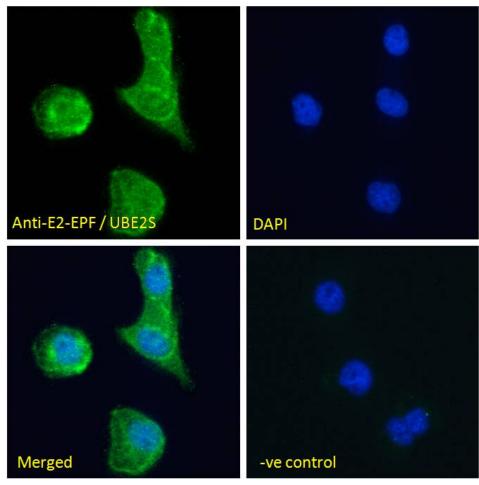

EB05758 Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic and nuclear staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

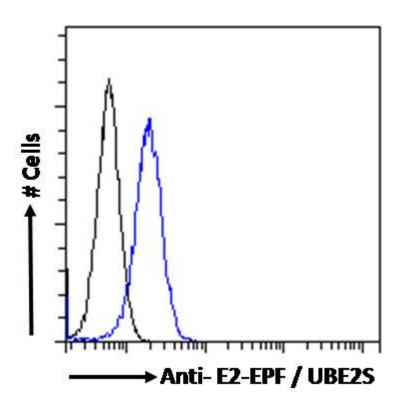

EB05758 Flow cytometric analysis of paraformaldehyde fixed HeLa cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.