

#### **International Office**

#### **Everest Biotech Ltd**

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

**Customer Service:** 

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB07120 - Goat Anti-Engrailed-2 Antibody

Size: 100µg specific antibody in 200µl



**Target Protein** 

Principal Names: Engrailed-2, engrailed homolog 2, EN2, AUTS1

Official Symbol: EN2

Accession Number(s): NP\_001418.2

Human GeneID(s): 2020

# Immunogen

Peptide with sequence C-NHSTTAKEGKSDSE, from the C Terminus of the protein sequence according to NP\_001418.2.

Please note the peptide is available for sale.

### **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

### **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:4000.


**Western blot:** Approx 35kDa band observed in Mouse Brain lysates (calculated MW of f 33.8kDa according to Mouse NP\_034264.1). Recommended concentration: 0.1-1µg/ml. Primary incubation 1 hour at room temperature.

**Immunofluorescence:** Strong expression of the protein seen in the nuclei of HEK293 and U2OS cells. Recommended concentration: 10µg/ml.


# **Species Reactivity**

Tested: Human, Mouse


Expected from sequence similarity: Human, Mouse, Rat, Dog

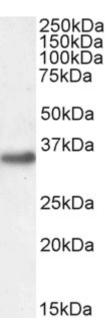

EB07120 (1μg/ml) staining of Mouse Brain lysate (35μg protein in RIPA buffer). Detected by chemiluminescence.

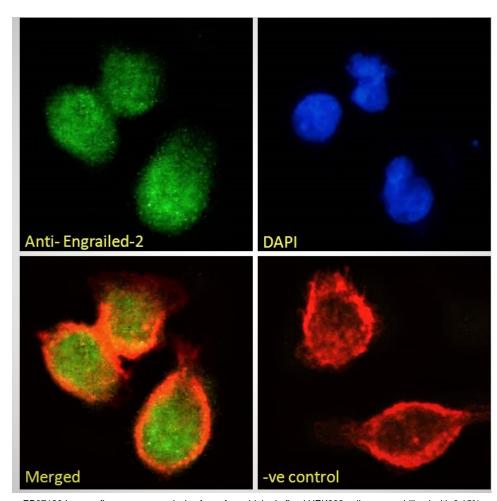

EB07120 Immunofluorescence analysis of paraformaldehyde fixed HEK293 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. Actin filaments were stained with phalloidin (red) and the nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

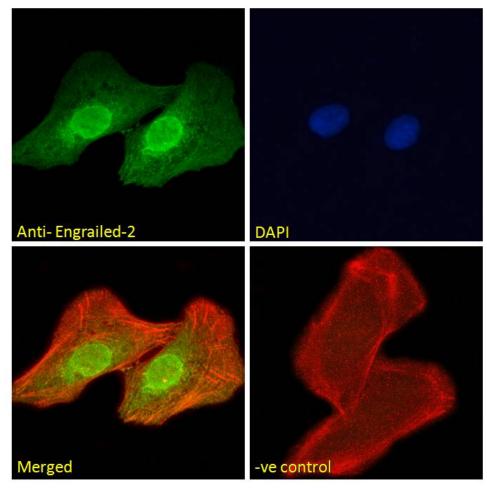

EB07120 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. Actin filaments were stained with phalloidin (red) and the nuclear stain is DAPI (blue). Negative control:

Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).