

#### **International Office**

#### **Everest Biotech Ltd**

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

**Customer Service:** 

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB06006-T - Goat Anti-FOXC1 (C Term) Antibody -

Size: 20µg specific antibody in 40µl

# **Target Protein**

**Principal Names:** FOXC1, forkhead box C1, ARA, IGDA, IHG1, FKHL7, IRID1, FREAC3, iridogoniodysgenesis type 1, forkhead (Drosophila)-like 7, forkhead-related activator 3, Forkhead, drosophila, homolog-like 7, FREAC-3, RIEG3, forkhead/winged helix-like

transcription factor 7, myeloid factor-delta

Official Symbol: FOXC1

Accession Number(s): NP\_001444.2

Human GeneID(s): 2296

Non-Human GeneID(s): 17300 (mouse)

### **Immunogen**

Peptide with sequence RTSGAFVYDCSKF, from the C Terminus of the protein sequence according to NP\_001444.2.

Please note the peptide is available for sale.

### **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

#### **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:32000.

Immunofluorescence: Strong expression of the protein seen in the cytoplasm and

nucleus of U2OS cells. Recommended concentration: 10µg/ml.


Flow Cytometry: Flow cytometric analysis of HEK293 cells. Recommended

concentration: 10ug/ml.


## **Species Reactivity**

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat

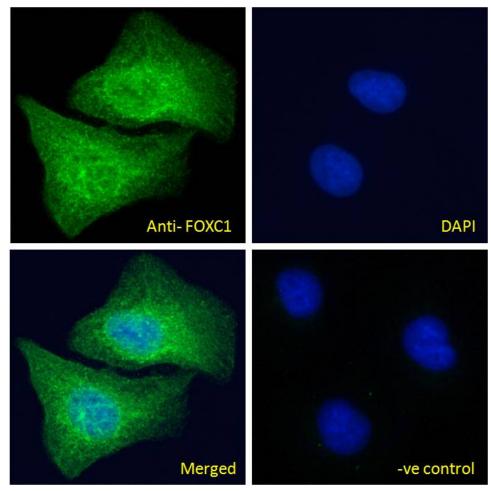

EB06006 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

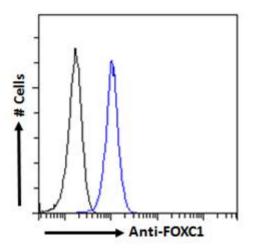

EB06006 Flow cytometric analysis of paraformaldehyde fixed HEK293 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.