



UK Office

Everest Biotech Ltd

Cherwell Innovation Centre
77 Heyford Park
Upper Heyford
Oxfordshire
OX25 5HD
UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326

www.everestbiotech.com

**Research Use Only. Not for
diagnostic or therapeutic use.**

EB08896 - Goat Anti-FURIN / PCSK3 Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: FURIN, furin (paired basic amino acid cleaving enzyme), FUR, PACE, PCSK3, SPC1, FES upstream region, dibasic processing enzyme, furin, furin, membrane associated receptor protein, paired basic amino acid cleaving enzyme (furin, membrane associated receptor protein), proprotein convertase subtilisin/kexin type 3

Official Symbol: FURIN

Accession Number(s): NP_002560.1

Human GeneID(s): [5045](#)

Immunogen

Peptide with sequence C-RTKRDVYQEPTDP, from the internal region of the protein sequence according to NP_002560.1.

Please note the [peptide](#) is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:32000.


Immunofluorescence: Strong expression of the protein seen in the plasma membrane and ER/cytoplasm of HeLa cells, and additionally the nucleus of U2OS cells.

Recommended concentration: 10µg/ml.


Species Reactivity

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat

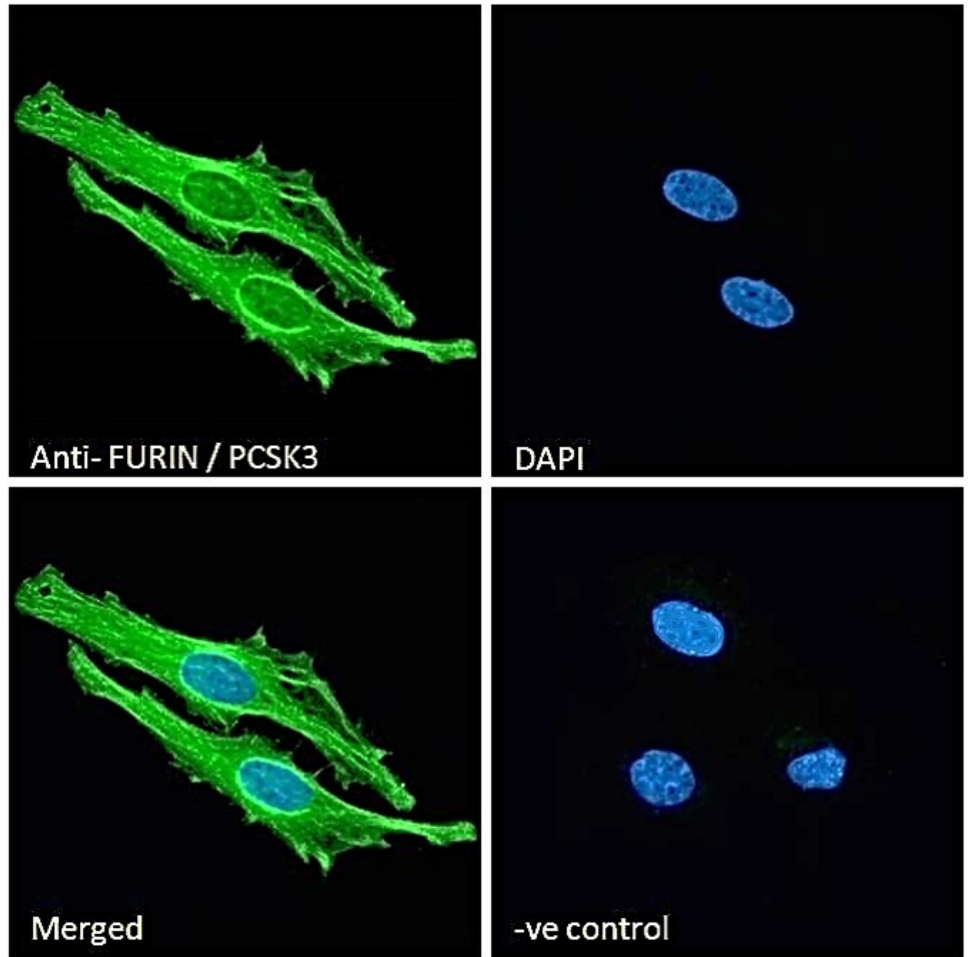

EB08896 Immunofluorescence analysis of paraformaldehyde fixed HeLa cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing plasma membrane and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

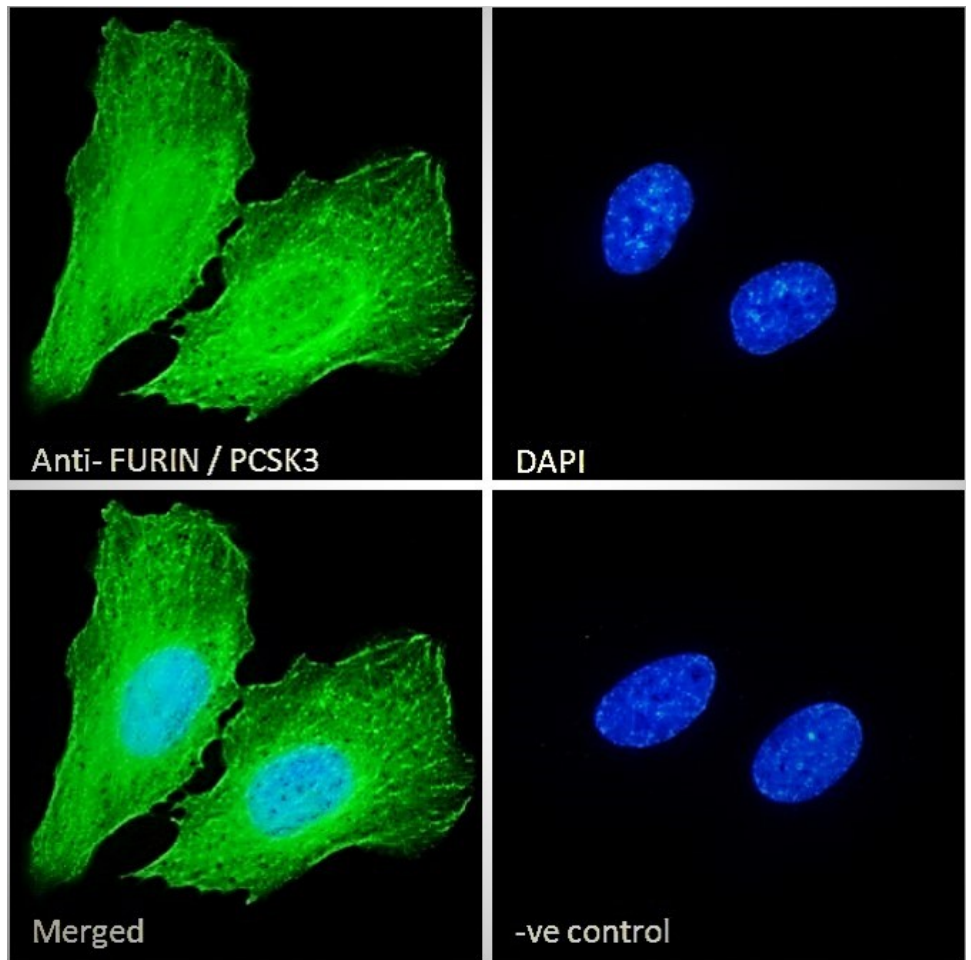

EB08896 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing plasma membrane, cytoplasmic and nuclear staining. The nuclear stain is DAPI (blue) Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).