



International Office

Everest Biotech Ltd

Vector Laboratories, Inc.
6737 Mowry Ave
Newark, CA 94560
United States

Customer Service:

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

**Research Use Only. Not for
diagnostic or therapeutic use.**

EB08514 - Goat Anti-IGFBP7 (aa194-207) Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: IGFBP7, insulin-like growth factor binding protein 7, FSTL2, IGFBP-7, IGFBP-7v, MAC25, PSF

Official Symbol: IGFBP7

Accession Number(s): NP_001544.1

Human GeneID(s): [3490](#)

Immunogen

Peptide with sequence C-NKVKRGHYGVQRTE, from the internal region of the protein sequence according to NP_001544.1 .

Please note the [peptide](#) is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:16000.


Western blot: Preliminary experiments gave an approx 50kDa band in Human Kidney, Spleen, Ovary, Prostate and Testis lysates after 0.5µg/ml antibody staining. Please note that currently we cannot find an explanation in the literature for the band we observe given the calculated size of 29.1kDa according to NP_001544.1. The 50kDa band was successfully blocked by incubation with the immunizing peptide.

IHC: Paraffin embedded Human Tonsil and Pancreas. Recommended concentration: 2.5µg/ml.


Species Reactivity

Tested: Human

Expected from sequence similarity: Human, Cow

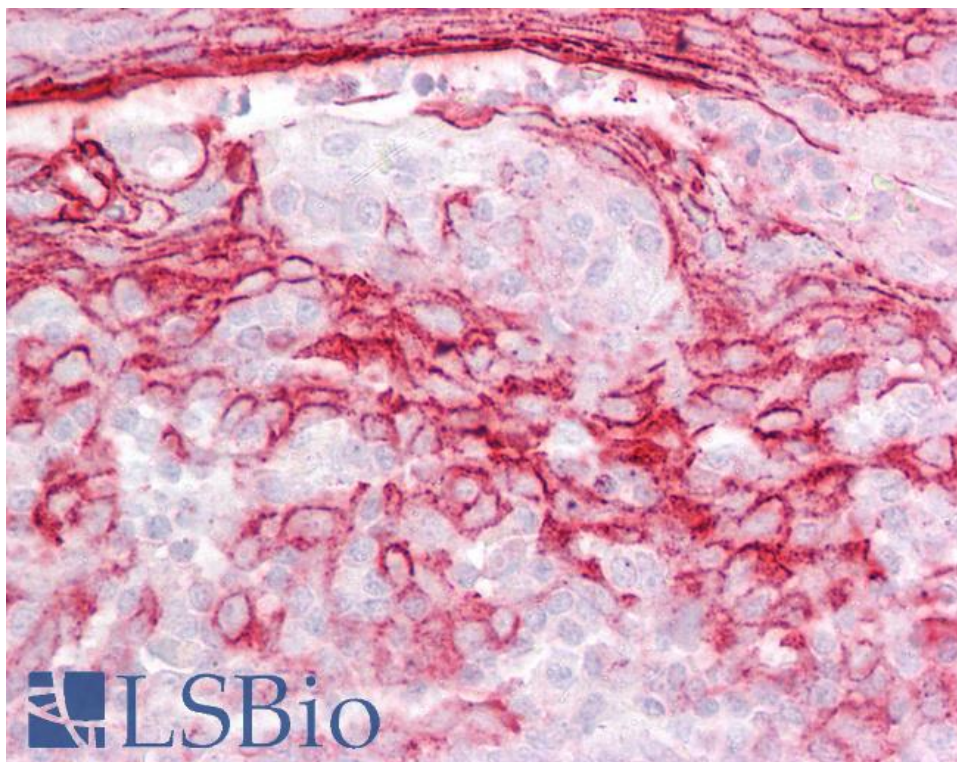

EB08514 (2.5µg/ml) staining of paraffin embedded Human Tonsil. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.

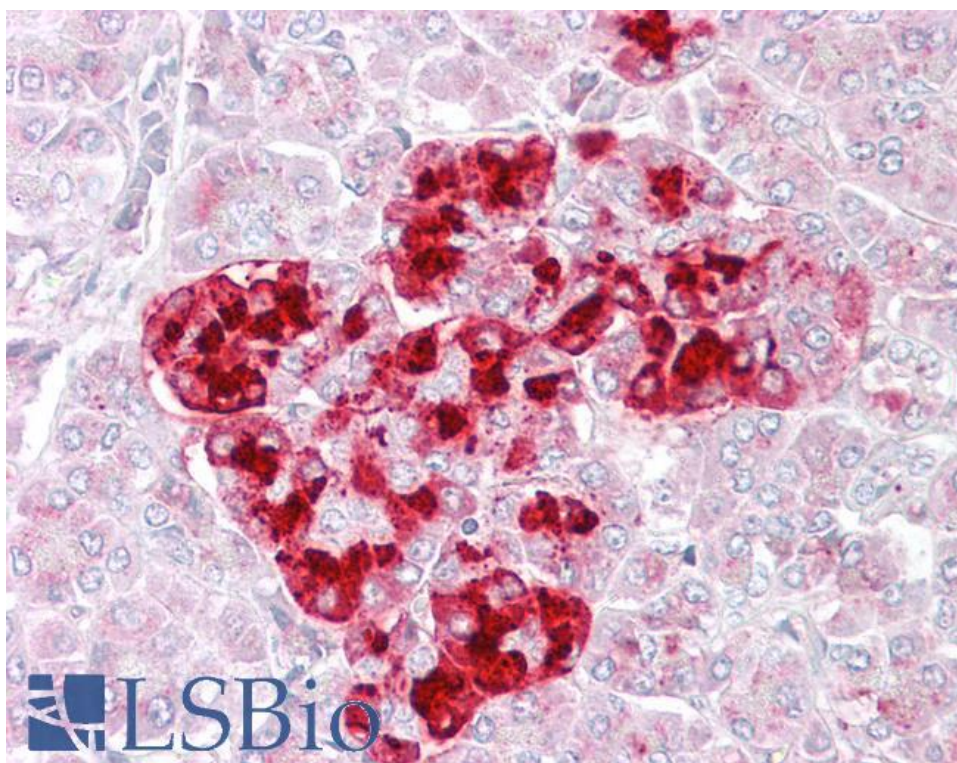

EB08514 (2.5µg/ml) staining of paraffin embedded Human Pancreas. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.