

International Office

Everest Biotech Ltd

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

Customer Service:

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB09583 - Goat Anti-IP6K3 Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: IP6K3, inositol hexakisphosphate kinase 3, IHPK3, INSP6K3, MGC102928, ATP:1D-myo-inositol-hexakisphosphate phosphotransferase, InsP6 kinase

3, OTTHUMP0000016212, inositol hexaphosphate kinase 3

Official Symbol: IP6K3

Accession Number(s): NP_473452.2

Human GeneID(s): 117283

Non-Human GenelD(s): 271424 (mouse), 688862 (rat)

Immunogen

Peptide with sequence RLCSEYPENKRHR, from the internal region of the protein sequence according to NP_473452.2.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:32000.


Immunofluorescence: Strong expression of the protein seen in HepG2 and U2OS cells. Recommended concentration: 10µg/ml.

Flow Cytometry: Flow cytometric analysis of HepG2 cells. Recommended concentration: 10ug/ml.


Species Reactivity

Tested: Human


Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

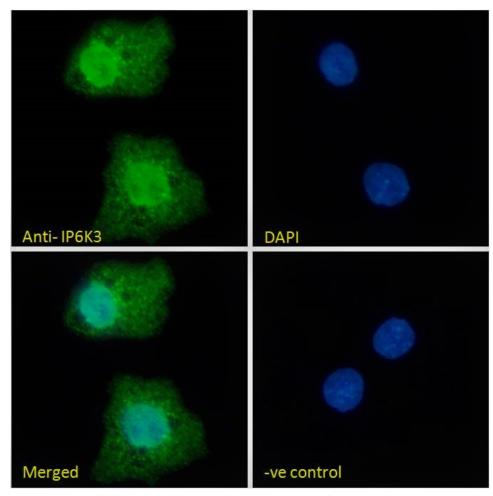

EB09583 Immunofluorescence analysis of paraformaldehyde fixed HepG2 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing strong nuclear and some ytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

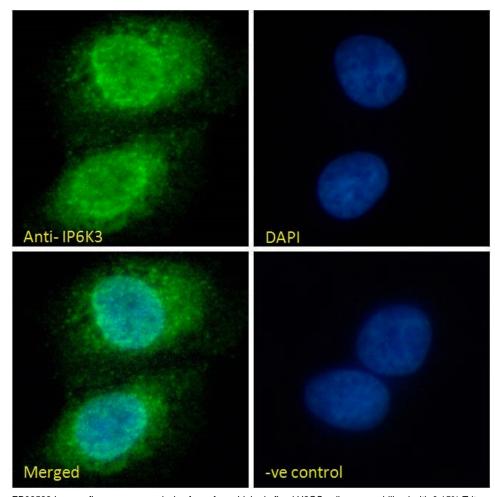

EB09583 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

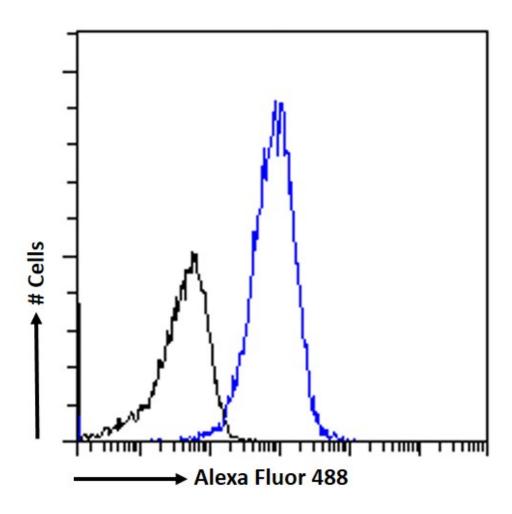

EB09583 Flow cytometric analysis of paraformaldehyde fixed HepG2 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.