

International Office

Everest Biotech Ltd

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

Customer Service:

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB05810-T - Goat Anti-KAP1 / TRIM28 Antibody - Trial

Size: 20µg specific antibody in 40µl



Target Protein

Principal Names: TRIM28, KAP1, TF1B, RNF96, TIF1B, tripartite motif-containing 28, KRAB-associated protein 1, nuclear corepressor KAP-1, transcriptional intermediary factor 1-beta, FLJ29029, KRAB-associated protein 1, tripartite motif-containing 28 protein

Official Symbol: TRIM28

Accession Number(s): NP_005753.1

Human GeneID(s): 10155

Immunogen

Peptide with sequence C-SSQELSGGPGDGP, from the C Terminus of the protein sequence according to NP_005753.1.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:2000.

IHC: In paraffin embedded Human Breast shows nuclear staining in epithelial duct cells. Recommended concentration: 2-4µg/ml.

Immunofluorescence: Strong expression of the protein seen in the nuclei of HeLa and U2OS cells. Recommended concentration: 10µg/ml.

Flow Cytometry: Flow cytometric analysis of HeLa cells. Recommended concentration: 10ug/ml.

Species Reactivity

Tested: Human

Expected from sequence similarity: Human, Mouse

Specific Reference


This antibody (previous batch) has been successfully used in IHC on Human:

Yokoe T, Toiyama Y, Okugawa Y, Tanaka K, Ohi M, Inoue Y, Mohri Y, Miki C, Kusunoki M


KAP1 Is Associated With Peritoneal Carcinomatosis in Gastric Cancer.

Ann Surg Oncol. 2010 Mar;17(3):821-8.


PMID: 19898899

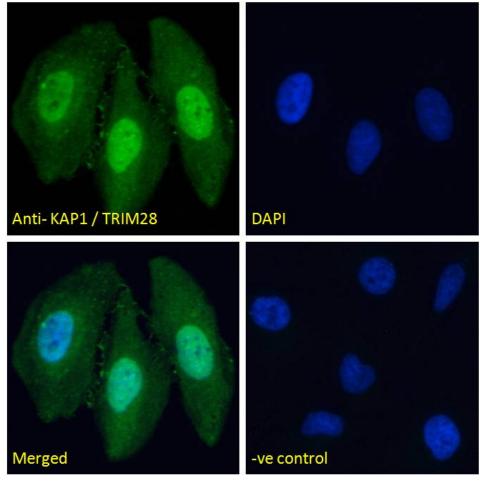

EB05810 Immunofluorescence analysis of paraformaldehyde fixed HeLa cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

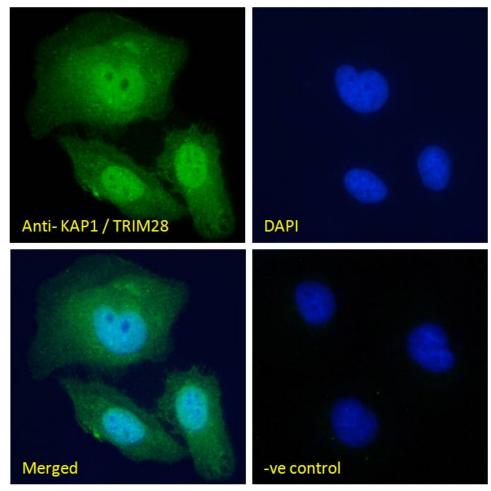


EB05810 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing nuclear staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

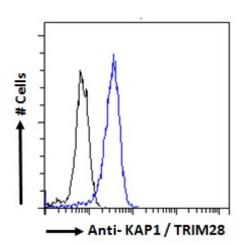

EB05810 Flow cytometric analysis of paraformaldehyde fixed HeLa cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.

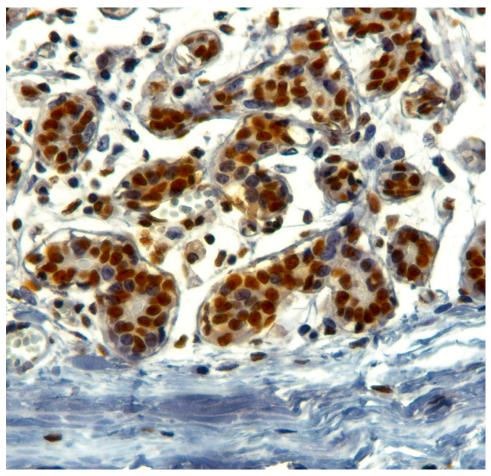

EB05810 ($2\mu g/ml$) staining of paraffin embedded Human Breast. Steamed antigen retrieval with Tris/EDTA buffer pH 9, HRP-staining.