

UK Office

Everest Biotech Ltd Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD UK

Enquiries: info@everestbiotech.com Sales: sales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

US Office

Everest Biotech c/o Abcore 405 Maple Street, Suite A106 Ramona, CA 92065 USA

Inquiries: info@everestbiotech.com Sales: usasales@everestbiotech.com Tech support: support@everestbiotech.com

Tel: 888-320-4628 (toll-free) Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB06664 - Goat Anti-Lipin 2 Antibody

Size: 100µg specific antibody in 200µl

Target Protein

Principal Names: LPIN2, KIAA0249, lipin 2 Official Symbol: LPIN2 Accession Number(s): NP_055461.1 Human GeneID(s): <u>9663</u>

Immunogen

Peptide with sequence CYWRDPIPEVDLDDLS, from the C Terminus of the protein sequence according to NP_055461.1.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:8000.


Immunofluorescence: Strong expression of the protein seen in the cytoplasm of A431 and U2OS cells. Recommended concentration: 10µg/ml.

Flow Cytometry: Flow cytometric analysis of A431 cells. Recommended concentration: 10ug/ml.


Species Reactivity

Tested: Human Expected from sequence similarity: Human, Mouse, Rat, Pig




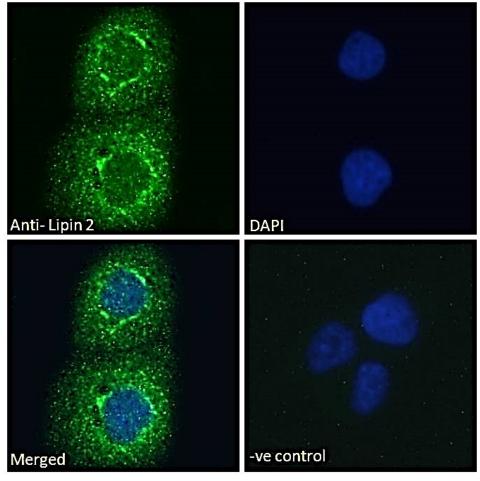

EB06664 Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

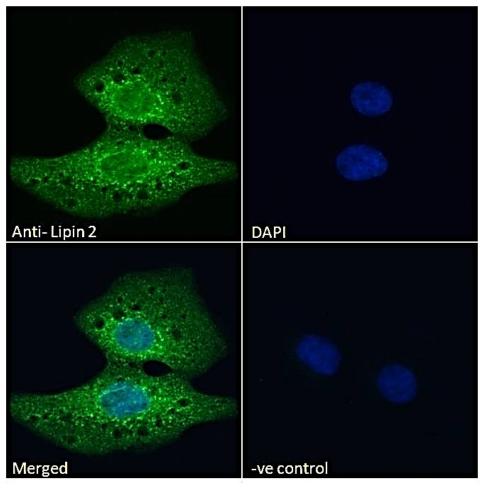

EB06664 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

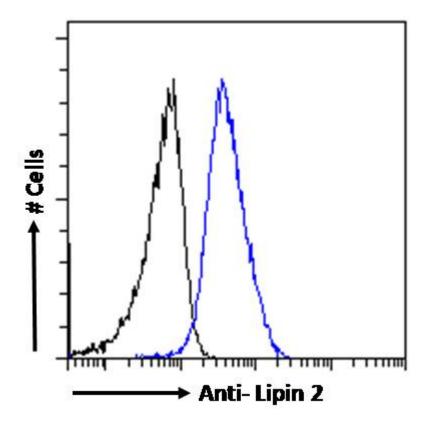

EB06664 Flow cytometric analysis of paraformaldehyde fixed A431 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control: Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.