

UK Office

Everest Biotech Ltd

Cherwell Innovation Centre

77 Heyford Park
Upper Heyford
Oxfordshire
OX25 5HD

UK

Enquiries:

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

US Office

Everest Biotech c/o Abcore

405 Maple Street, Suite A106

Ramona, CA 92065 USA

Inquiries:

info@everestbiotech.com

Sales:

 $\underline{usasales@everest biotech.com}$

Tech support:

support@everestbiotech.com

Tel: 888-320-4628 (toll-free)

Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB06289 - Goat Anti-MDG1 / DNAJB9 Antibody

Size: 100µg specific antibody in 200µl



Target Protein

Principal Names: DNAJB9, MDG1, ERdj4, DKFZP564F1862, DnaJ (Hsp40) homolog, subfamily B, member 9, endoplasmic reticulum DnaJ homolog 4, microvascular

endothelial differentiation gene 1, MST049, MSTP049

Official Symbol: DNAJB9

Accession Number(s): NP_036460.1

Human GeneID(s): 4189

Non-Human GenelD(s): 27362 (mouse), 24908 (rat)

Immunogen

Peptide with sequence RGNMVTTYTDCSGQ, from the C Terminus of the protein sequence according to NP_036460.1.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested

Peptide ELISA: antibody detection limit dilution 1:4000.

Western blot: Preliminary experiments gave no signal but low background in A459 and

Human Skeletal Muscle extracts at up to 1µg/ml.


IHC: In paraffin embedded Human Kidney shows ER staining in PCT. Recommended

concentration: 4-6µg/ml.

Species Reactivity

Tested: Human

Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

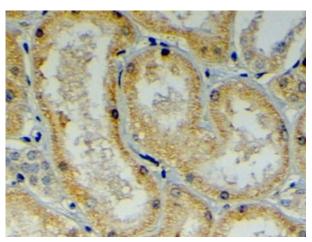

EB06289 (4µg/ml) staining of paraffin embedded Human Kidney. Steamed antigen retrieval with citrate buffer pH 6, HRP-staining.