



## UK Office

### Everest Biotech Ltd

Cherwell Innovation Centre  
77 Heyford Park  
Upper Heyford  
Oxfordshire  
OX25 5HD  
UK

Enquiries:

[info@everestbiotech.com](mailto:info@everestbiotech.com)

Sales:

[sales@everestbiotech.com](mailto:sales@everestbiotech.com)

Tech support:

[support@everestbiotech.com](mailto:support@everestbiotech.com)

Tel: +44 (0)1869 238326

[www.everestbiotech.com](http://www.everestbiotech.com)

**Research Use Only. Not for  
diagnostic or therapeutic use.**

## EB07007 - Goat Anti-Netrin 1 Antibody

Size: 100µg specific antibody in 200µl



### Target Protein

**Principal Names:** netrin 1, HGNC:8029, NTN1, NTN1L, Netrin-1, netrin 1, mouse, homolog of

**Official Symbol:** NTN1

**Accession Number(s):** NP\_004813.2

**Human GeneID(s):** [9423](#)

**Non-Human GeneID(s):** 18208 (mouse), 114523 (rat)

### Immunogen

Peptide with sequence KFQQREKKGKCKKA, from the C Terminus of the protein sequence according to NP\_004813.2.

Please note the [peptide](#) is available for sale.

### Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

### Applications Tested

**Peptide ELISA:** antibody detection limit dilution 1:32000.

**Immunofluorescence:** Strong expression of the protein seen in the membrane and cytoplasm of U2OS cells. Recommended concentration: 10µg/ml.

### Species Reactivity

**Tested:** Human

**Expected from sequence similarity:** Human, Mouse, Rat, Dog, Cow

### Specific References

**This antibody (previous batch) has been successfully used in Western blot on Human:**

Jayakumar C, Nauta FL, Bakker SJ, Bilo H, Gansevoort RT, Johnson MH, Ramesh G. Netrin-1, a urinary proximal tubular injury marker, is elevated early in the time course of human diabetes.

J Nephrol. 2014 Apr;27(2):151-7.

PMID: 24510764


**This antibody (previous batch) has been successfully used in Western blot on Human and Mouse:**

White JJ, Mohamed R, Jayakumar C, Ramesh G.

Tubular injury marker netrin-1 is elevated early in experimental diabetes.

J Nephrol. 2013 Nov-Dec;26(6):1055-64.

PMID: 24052471

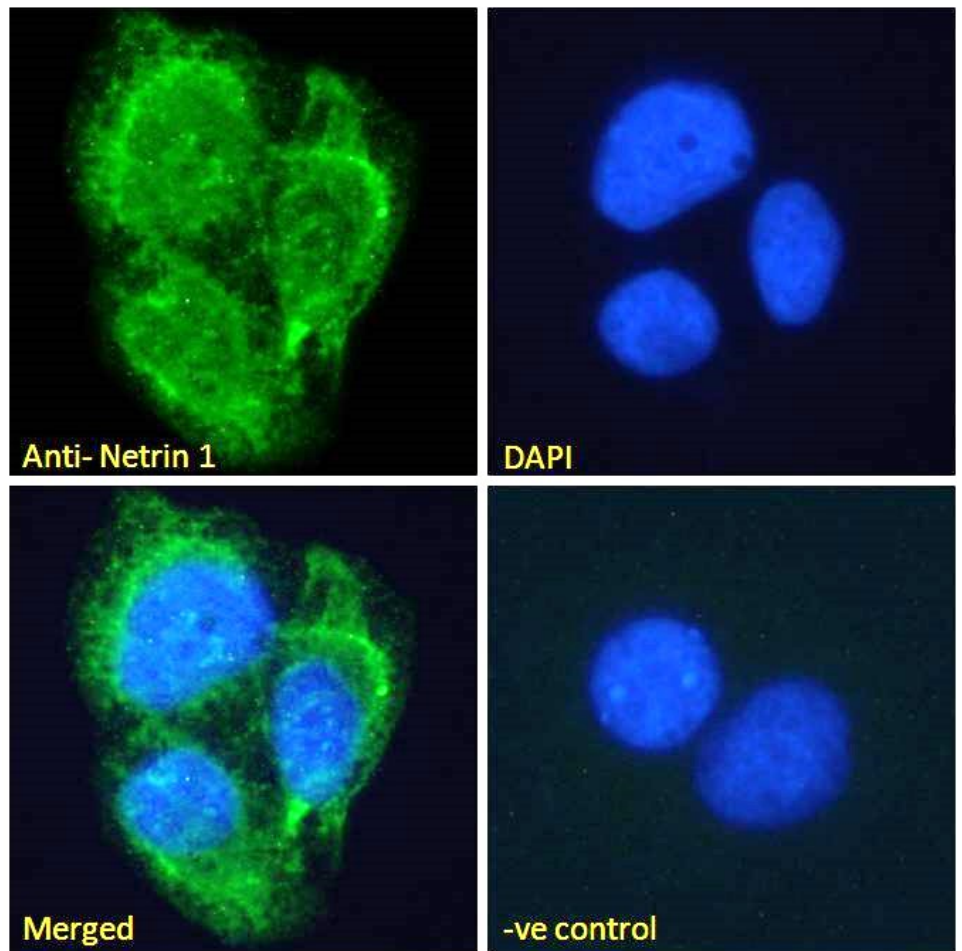

EB07007 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing membrane and secreted staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).