

International Office

Everest Biotech Ltd Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

Customer Service: <u>customerservice@vectorlabs.com</u> Technical Service: <u>technical@vectorlabs.com</u>

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

EB06660 - Goat Anti-Plexin D1 Antibody

Size: 100µg specific antibody in 200µl

C

Target Protein Principal Names: PLXND1, PLEXD1, KIAA0620, plexin D1, MGC75353 Official Symbol: PLXND1 Accession Number(s): NP_055918.2 Human GenelD(s): 23129

Immunogen

Peptide with sequence C-LAEPKKSHRQSH, from the internal region of the protein sequence according to NP_055918.2.

Please note the peptide is available for sale.

Purification and Storage

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

Applications Tested


Peptide ELISA: antibody detection limit dilution 1:16000. IHC: .

Immunofluorescence: Strong expression of the protein seen in the plasma membrane and cytoplasm of HeLa and U251cells. Recommended concentration: 10µg/ml.


Flow Cytometry: Flow cytometric analysis of K562 cells. Recommended concentration: 10ug/ml.

Species Reactivity


Tested: Human Expected from sequence similarity: Human, Dog, Pig, Cow

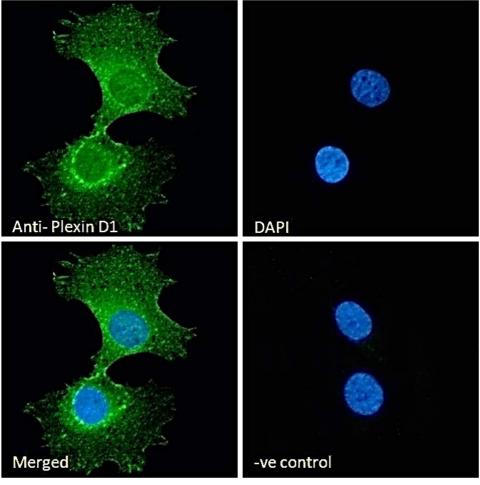

EB06660 Immunofluorescence analysis of paraformaldehyde fixed HeLa cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing membrane and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

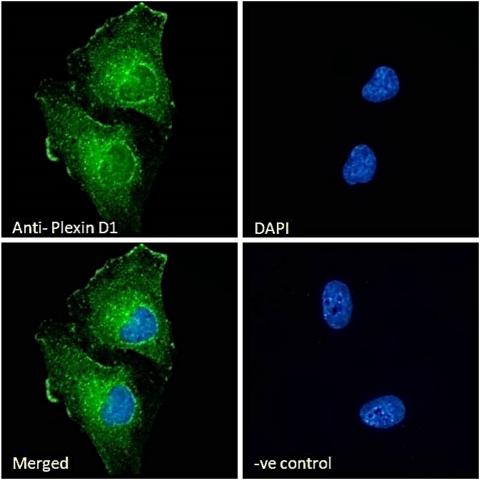

EB06660 Immunofluorescence analysis of paraformaldehyde fixed U251 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing membrane and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

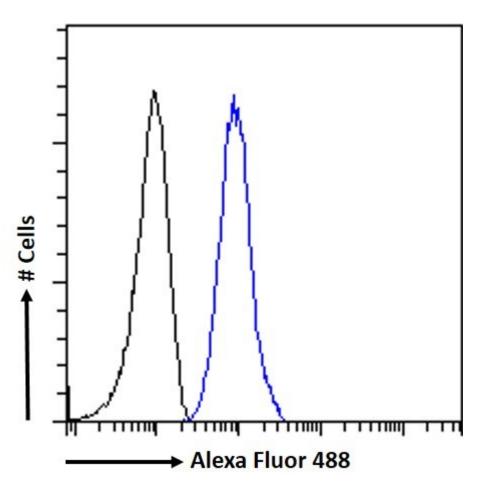

EB06660 Flow cytometric analysis of paraformaldehyde fixed K562 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control: Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.