

## **UK Office**

#### **Everest Biotech Ltd**

Cherwell Innovation Centre 77 Heyford Park Upper Heyford Oxfordshire OX25 5HD

UK

**Enquiries:** 

info@everestbiotech.com

Sales:

sales@everestbiotech.com

Tech support:

support@everestbiotech.com

Tel: +44 (0)1869 238326 Fax: +44 (0)1869 238327

#### **US Office**

#### **Everest Biotech c/o Abcore**

405 Maple Street, Suite A106

Ramona, CA 92065 USA

Inquiries:

info@everestbiotech.com

Sales:

 $\underline{usasales@everest biotech.com}$ 

Tech support:

support@everestbiotech.com

Tel: 888-320-4628 (toll-free)

Fax: 888-841-9041

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB06660 - Goat Anti-Plexin D1 Antibody

Size: 100µg specific antibody in 200µl



## **Target Protein**

Principal Names: PLXND1, PLEXD1, KIAA0620, plexin D1, MGC75353

Official Symbol: PLXND1

Accession Number(s): NP\_055918.2

Human GenelD(s): 23129

#### **Immunogen**

Peptide with sequence C-LAEPKKSHRQSH, from the internal region of the protein sequence according to NP\_055918.2.

Please note the peptide is available for sale.

## **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

## **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:16000.

IHC: .


**Immunofluorescence:** Strong expression of the protein seen in the plasma membrane and cytoplasm of HeLa and U251cells. Recommended concentration: 10µg/ml.

**Flow Cytometry:** Flow cytometric analysis of K562 cells. Recommended concentration: 10ug/ml.


# **Species Reactivity**

Tested: Human


Expected from sequence similarity: Human, Dog, Pig, Cow

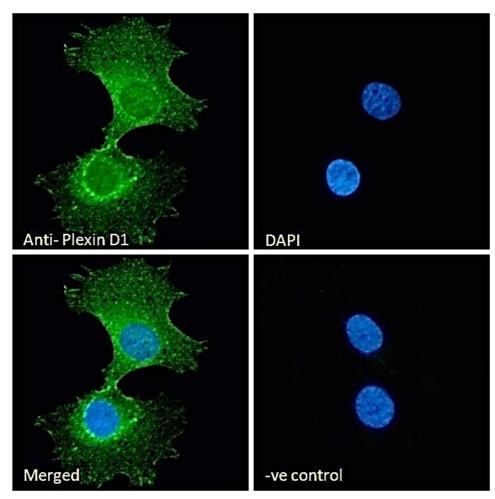

EB06660 Immunofluorescence analysis of paraformaldehyde fixed HeLa cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing membrane and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

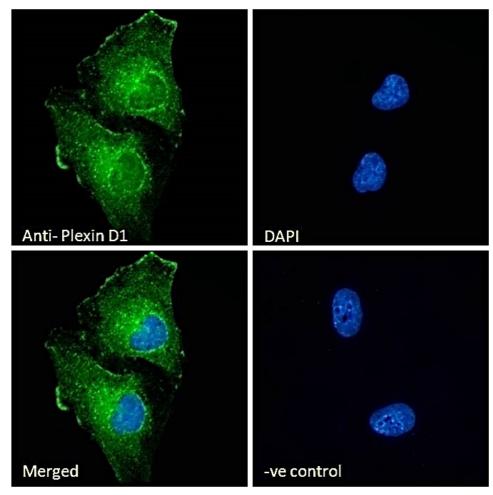

EB06660 Immunofluorescence analysis of paraformaldehyde fixed U251 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing membrane and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

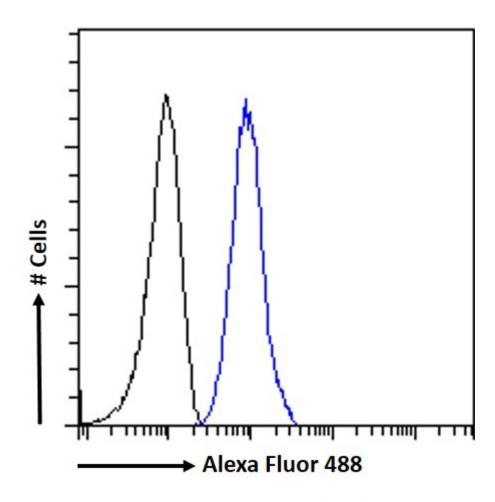

EB06660 Flow cytometric analysis of paraformaldehyde fixed K562 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.