

#### International Office

#### **Everest Biotech Ltd**

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

**Customer Service:** 

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB09512 - Goat Anti-Restin / CLIP1 Antibody

Size: 100µg specific antibody in 200µl



### **Target Protein**

**Principal Names:** CLIP1, CAP-GLY domain containing linker protein 1, CLIP, CLIP-170, CLIP170, CYLN1, MGC131604, RSN, Reed-Steinberg cell-expressed intermediate filament-associated protein, cytoplasmic linker 1, restin, restin (Reed-Steinberg cell-expressed intermediate filament-associated protein)

Official Symbol: CLIP1

Accession Number(s): NP\_002947.1; NP\_937883.1; NP\_001234926.1;

NP\_001376220.1

Human GeneID(s): 6249

Non-Human GenelD(s): 56430 (mouse), 65201 (rat)

Important Comments: This antibody is expected to recognize both reported isoforms

(NP\_002947.1; NP\_937883.1).

#### **Immunogen**

Peptide with sequence C-PHSTHHGSRGEERP, from the C Terminus of the protein sequence according to NP\_002947.1; NP\_937883.1; NP\_001234926.1; NP\_001376220.1.

Please note the peptide is available for sale.

### **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin.

Aliquot and store at -20°C. Minimize freezing and thawing.

## **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:8000.

**Western blot:** Approx. 150kDa band observed in lysates of cell lines HEK293, A549, HepG2 and HeLa (calculated MW of 161kDa according to NP\_002947.1). Recommended concentration: 0.1-1µg/ml. Primary incubation 1 hour at room temperature.

IHC: Paraffin embedded Human Thyroid and Prostate. Recommended concentration: 2.5µg/ml.


**Immunofluorescence:** Strong expression of the protein seen in A431 and A549 cells. Recommended concentration: 10µg/ml.

**Flow Cytometry:** Flow cytometric analysis of A549 cells. Recommended concentration: 10ug/ml.

### **Species Reactivity**

Tested: Human


Expected from sequence similarity: Human

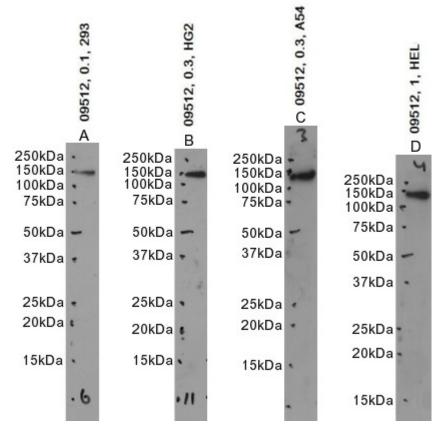


EB09512 optimised QC. Primary incubation 1 hour at room temperature.

Image A: HEK293 cell lysate at primary Ab concentration 0.1ug/ml, Images B, C: HepG2, A549 cell lysate at primary Ab concentration 0.3ug/ml, and Image D: HeLa cell lysate at primary Ab concentration 1ug/ml. (Loaded 35µg protein in RIPA buffer, per lane).


Detected by chemiluminescence.

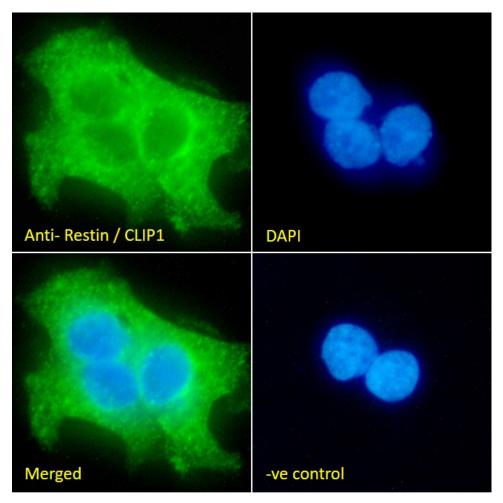

EB09512 Immunofluorescence analysis of paraformaldehyde fixed A431 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

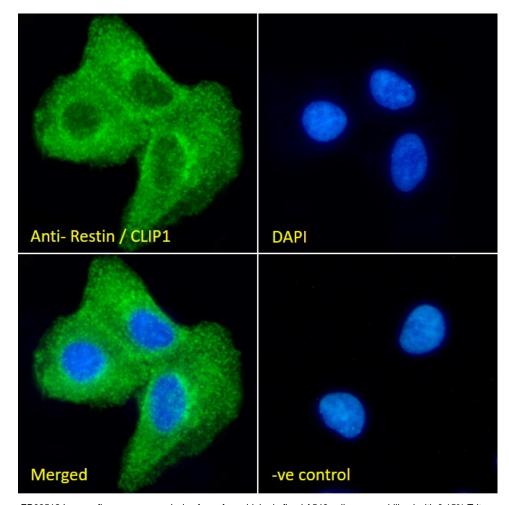


EB09512 Immunofluorescence analysis of paraformaldehyde fixed A549 cells, permeabilized with 0.15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml), showing cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

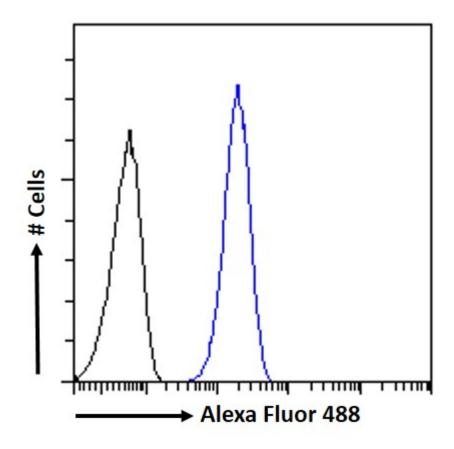


EB09512 Flow cytometric analysis of paraformaldehyde fixed A549 cells (blue line), permeabilized with 0.5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control:

Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.

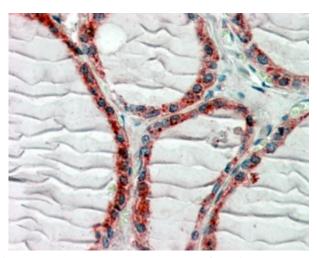

 $EB09512\ (2.5\mu g/ml)\ staining\ of\ paraffin\ embedded\ Human\ Thyroid\ Gland.\ Steamed\ antigen\ retrieval\ with\ citrate$   $buffer\ pH\ 6,\ AP-staining.$ 

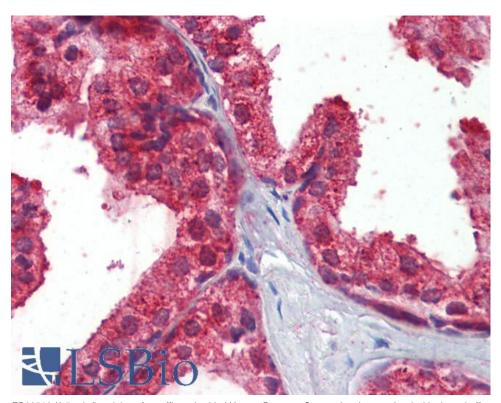

EB09512 (2.5μg/ml) staining of paraffin embedded Human Prostate. Steamed antigen retrieval with citrate buffer Ph 6, AP-staining.