

### **International Office**

#### **Everest Biotech Ltd**

Vector Laboratories, Inc. 6737 Mowry Ave Newark, CA 94560 United States

**Customer Service:** 

customerservice@vectorlabs.com

Technical Service:

technical@vectorlabs.com

Tel: +1 (800) 227-6666

www.everestbiotech.com

Research Use Only. Not for diagnostic or therapeutic use.

# EB06831 - Goat Anti-VPS41 Antibody

Size: 100µg specific antibody in 200µl



### **Target Protein**

**Principal Names:** VPS41, vacuolar protein sorting 41 (yeast), HGNC:12713, HVSP41, vacuolar assembly protein 41, vacuolar protein sorting 41 (yeast homolog), vacuolar protein sorting 41 homolog (S. cerevisiae), HVPS41, hVps41p, vacuolar protein sorting 41

Official Symbol: VPS41

Accession Number(s): NP\_055211.2; NP\_542198.2

Human GeneID(s): 27072

Important Comments: This antibody is expected to recognise both reported isoforms

(NP\_055211.2 and NP\_542198.2)

### **Immunogen**

Peptide with sequence C-DHLKKDSQNKTLLK, from the internal region of the protein sequence according to NP\_055211.2; NP\_542198.2.

Please note the peptide is available for sale.

# **Purification and Storage**

Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.

Supplied at 0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin

Aliquot and store at -20°C. Minimize freezing and thawing.

# **Applications Tested**

Peptide ELISA: antibody detection limit dilution 1:16000.


**Western blot:** Preliminary experiments gave bands at approx 48kDa and 28kDa in human duodenum, kidney and placenta lysates after 1µg/ml antibody staining. Please note that currently we cannot find an explanation in the literature for the bands we observe given the calculated size of 98.5kDa according to NP\_055211.1 and 100kDa according to NP\_542198.1. Both detected bands were successfully blocked by incubation with the immunizing peptide (and BLAST results with the immunizing peptide sequence did not identify any other proteins to explain the additional bands).

**IHC:** In paraffin embedded Human Kidney shows staining of collective tubules. Recommended concentration: 4-6μg/ml. Paraffin embedded Human Placenta and Prostate. Recommended concentration: 3.75μg/ml.


#### **Species Reactivity**

Tested: Human


Expected from sequence similarity: Human, Mouse, Rat, Dog, Cow

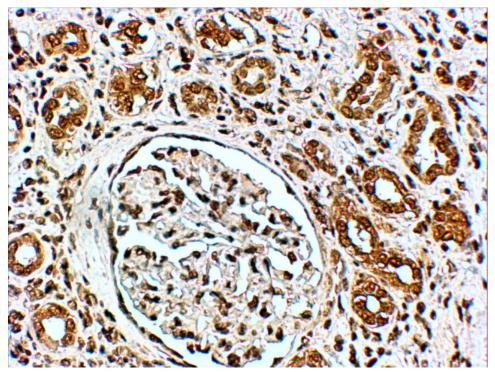

EB06831 (4µg/ml) staining of paraffin embedded Human Kidney. Steamed antigen retrieval with citrate buffer pH 6, HRP-staining. Similar results were obtained after antigen retrieval at pH9.

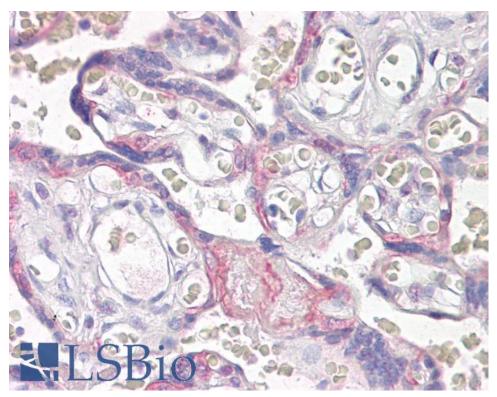

EB06831 (3.75µg/ml) staining of paraffin embedded Human Placenta. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.

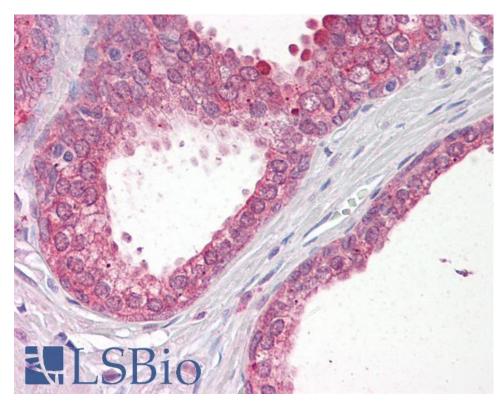

EB06831 (3.75 $\mu$ g/ml) staining of paraffin embedded Human Prostate. Steamed antigen retrieval with citrate buffer pH 6, AP-staining.