



## UK Office

### Everest Biotech Ltd

Cherwell Innovation Centre  
77 Heyford Park  
Upper Heyford  
Oxfordshire  
OX25 5HD  
UK

Enquiries:

[info@everestbiotech.com](mailto:info@everestbiotech.com)

Sales:

[sales@everestbiotech.com](mailto:sales@everestbiotech.com)

Tech support:

[support@everestbiotech.com](mailto:support@everestbiotech.com)

Tel: +44 (0)1869 238326

[www.everestbiotech.com](http://www.everestbiotech.com)

**Research Use Only. Not for  
diagnostic or therapeutic use.**

## SB-SHIFIX50PLUS - SB-SHIFIX50PLUS Coverslips

Size: Pack of 50


### Shi-fix™ coverslips

Do you wish immunofluorescence of suspension cells was as easy as with adherent cells? Now it is! Shi-fix™ (manufactured by Shikhar Biotech) cover-slips allow suspension cells to be layered on them or directly grown as a monolayer. Easy, hassle-free protocol and no spin columns needed for fixing suspension cells to slides. Just add your cells to the coverslips, wait for 30 mins, wash unbound cells with PBS and proceed to immunostaining (or continue the culture to obtain desired cell densities).


These cover-slips can also be used for staining suspension cells for microscopy, or for fixing suspension cells for confocal microscopy. What's more, you can also stimulate the cells while they are attached on the cover-slip, should you need to do so. You can work on non-adherent cells as if they were adherent!

### Suspension cell immunofluorescence with Shi-fix™ coverslips


**SIZE:** ø 15mm x 0.17mm(#1.5)

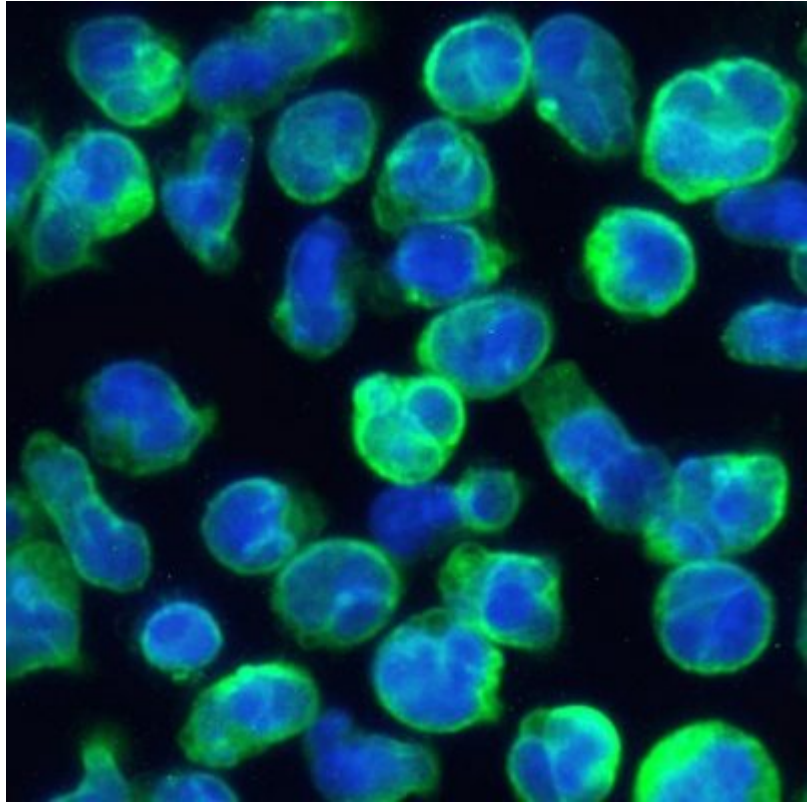


Immunofluorescence analysis of Molt4 cells immobilized on SHI-FIX™ coverslips.

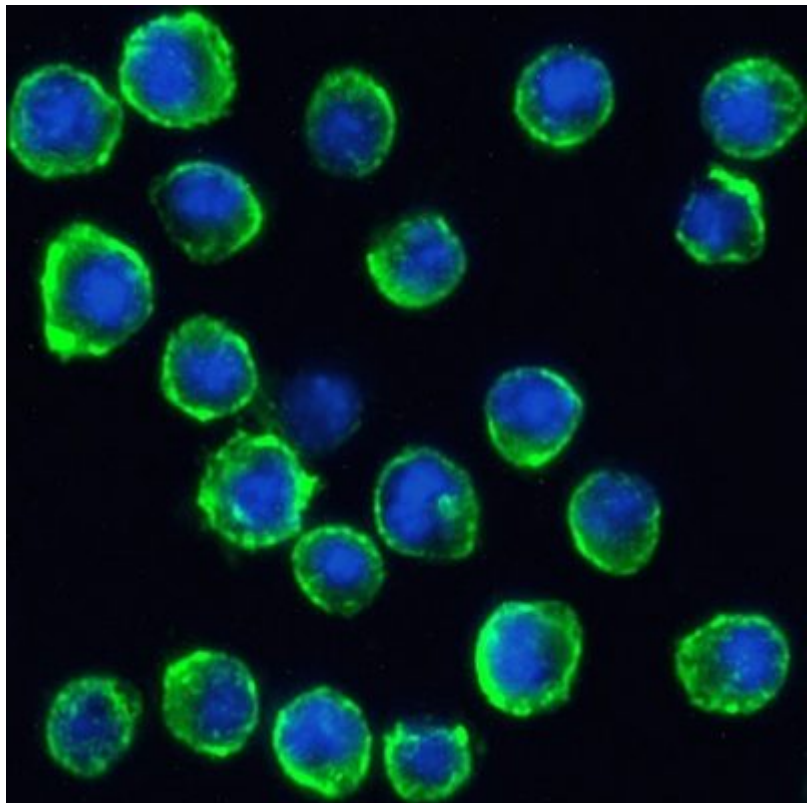

Immunofluorescence analysis of Daudi cells immobilized on SHI-FIX™ coverslips.

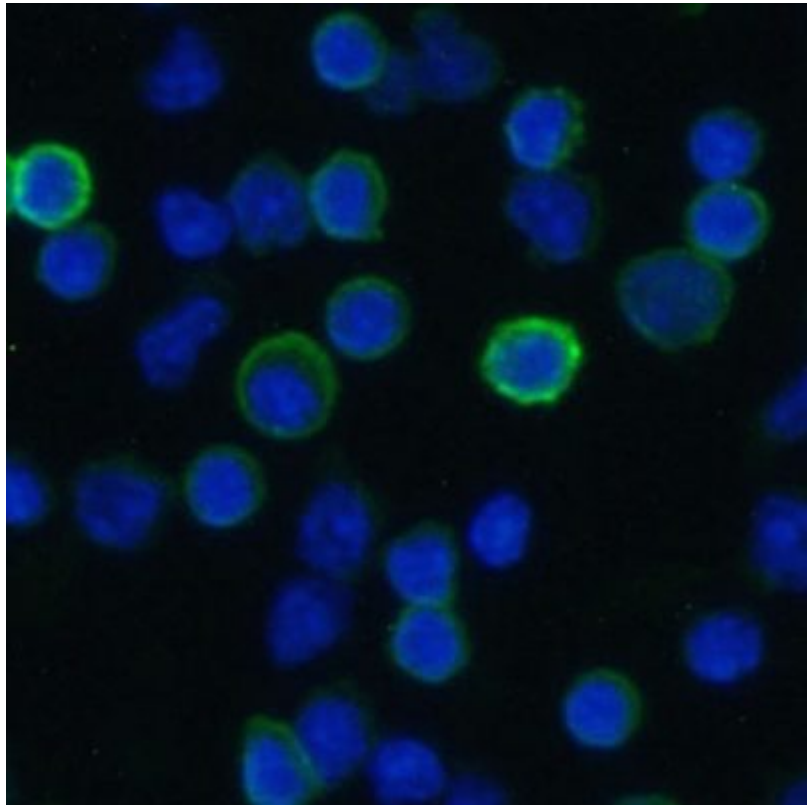

Immunofluorescence analysis of rat splenocytes immobilized on SHI-FIX™ coverslips

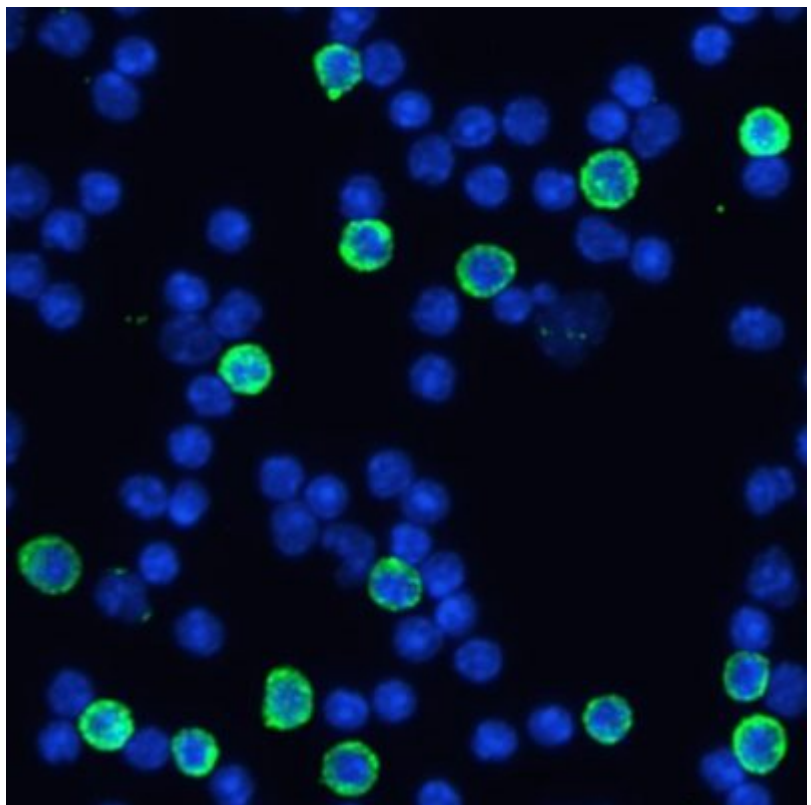

Immunofluorescence analysis of mouse thymocytes immobilized on SHI-FIX™ coverslips